EM-100

## EN User Guide

Thank you for purchasing this Lifetrons® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

### A. Insert the Battery

1. Turn the device to the back side. Use your fingers to press and slide the back cover downward. Then lift and remove the back cover.
2. Insert 1 AAA battery into the battery compartment. Be sure to observe the polarity as illustrated in the drawing ①.
3. Slide the cover back on to close the battery compartment and ensure that the closure latch has been securely closed.

### B. How to Use

1. Always use on clean skin. Apply a thin layer of serum, moisturizer or cream to your eye area.
2. Hold the device with your hand touching the conductive pole strip to activate the negative ion energy. Gently press the massage head onto your skin to start the eye massage. **Slowly and gently** move the device over all points around your eyes with light pressure.
3. For best results, use up to 4 times a week for 10 minutes per session.
4. The device will automatically stop massaging when pressure is no longer applied. After use, gently clean the massage head with a cotton pad.

**Note:** Skin must be clean before use. Serum, moisturizer or cream (not included) must be applied in combination with this device to ensure proper operation and maximum benefits.

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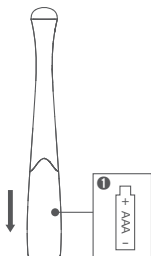
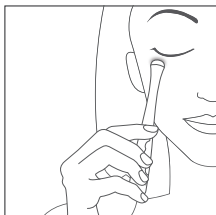
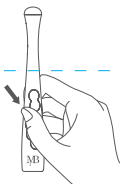
# EN Safety Precautions

## Safety Precautions

Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

- Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
- Do not modify, repair or dismantle the product. Doing so may result in fire, electric shock, complete breakdown of the product, etc. and will void the warranty.
- Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
- Do not use this product in locations with high temperatures or high humidity levels (e.g., bathrooms); ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
- Do not insert foreign objects or liquids into the product; this could damage the product.
- Refrain from dropping the product or bumping it against other objects.
- Do not use this product in hazardous locations.
- Should any abnormal situations or problems occur while using this product, discontinue its use immediately.
- The product may explode if disposed of in a fire. Immediately discontinue using this device if an allergic reaction occurs or you feel discomfort.
- This product is not suitable to use on body parts that are embedded with silicone, metal or other foreign objects. Do not use directly on the eyes.
- It is recommended to wait at least 6 months after any facial treatments such as a face lift or Botox injections before using this product. If you are unsure, consult with your doctor.
- This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.
- Do not place this product near strong magnetic or electric fields - this may cause product failure to occur.
- This product has a special lacquer - when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.

## Children

Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

## Disposal of Electrical and Electronic Equipment

This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, household waste disposal service or where you purchased the product.

## Limited Warranty

If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

1. The original invoice or sales receipt with the date of purchase is presented with the defective unit.
2. Lifetrons does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
3. The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrons authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrons reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
4. The warranty does not cover product failures which have been caused by use of accessories or other peripheral devices which are not Lifetrons brand original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental breakage.
5. The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
6. Any defective unit or part shall become property of the manufacturer.

## Additional Information

Please visit our website [www.lifetrons.com](http://www.lifetrons.com) or contact our customer service center at [cs@lifetrons.com](mailto:cs@lifetrons.com) to learn more about our new products and special offers.

Product specifications are subject to change at the sole decision and discretion of Lifetrons Switzerland AG and Lifetrons International Limited without prior notice. Any devices or accessories not listed under "What's Inside" are in reference to compatibility only and are not included. Product images are for reference purposes only, color may vary from the actual product.

"Lifetrons" is a trademark or registered trademark owned by Lifetrons Switzerland AG or Lifetrons International Limited. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrons Switzerland AG or Lifetrons International Limited. All rights are reserved.

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

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