


## Ionic Facial & Body Spa


4-in-1 with Microcurrent Technology

EP-500B

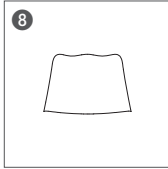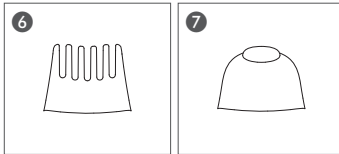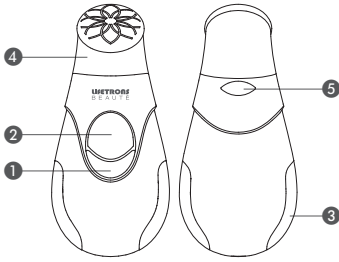C

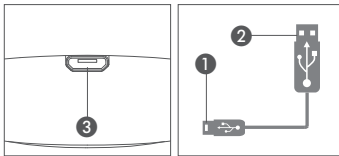D

## EN User Guide

### User Guide

Thank you for purchasing this Lifetrons® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

### How to Use

#### A. Changing the Massage Head

1. Press the release button **1** on the back of the device to release the massage head.
2. Select the massage head you would like to use and attach it **4** / **6** / **7** / **8** to the body of the device.

**Note:** Do not change the massage head when the power is on.

#### B. Treatment Methods

1. Start by attaching the facial massage head **1** to the device. Press the power button **1** 1 time to turn on the device. The default program will begin with the facial massage in negative ion mode. Hold the device with both sides of your hand touching the conductive pole strips **2** to activate the negative ion energy while gently massaging the face. The device will automatically shut down after 2 minutes of continuous use.
2. Press **1** 2 times to begin the positive ion cleansing mode with the facial massage head **1**. Hold the device with your hand touching the conductive pole strips **2** to activate the positive ion energy while gently massaging the face. The device will automatically shut down after 3 minutes of continuous use.
3. Attach the eye massage head **3** to the device. Press **1** 3 times for the eye massage in negative ion mode. Hold the device with your hand touching the conductive pole strips **2** to activate the negative ion energy while gently massaging around the eyes. The device will automatically shut down after 3 minutes of continuous use.
4. Attach the body massage head **2** to the device. Press **1** 4 times for the body massage in negative ion mode. Hold the device with your hand touching the conductive pole strips **2** to activate the negative ion energy while gently massaging the body. The device will automatically shut down after 5 minutes of continuous use.
5. Attach the scalp massage head **8** to the device. Press **1** 5 times for the scalp massage in negative ion mode. Hold the device with your hand touching the conductive pole strips **2** to activate the negative ion energy while gently gliding the massager over your entire scalp. The device will automatically shut down after 2 minutes of continuous use.

**Note:** Use on clean skin and be sure to remove makeup before use. This product is equipped with an automatic detection function. When the massage head and the conductive pole strips are not connected to the human body, the device will beep 3 times every 5 seconds and completely shut down after 15 seconds. If you want to shorten the usage time, simply press and hold the power button for 3 seconds. The device will beep 3 times before shutting down.

#### C. Recharging the Device

1. Connect part **1** (micro-USB jack) to part **3** (micro-USB port).
2. Connect part **2** (USB jack) to a USB port with power supply (e.g., AC adaptor with USB power supply).

#### D. Battery Indicator

**Low battery:** The LCD displays a flashing battery icon D (only when massager is in use).  
**When charging:** The LCD displays a rolling battery icon.  
**Fully charged battery:** The LCD displays nothing.

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## EN Safety Precautions

### Safety Precautions

Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

- Discontinue using the product immediately if any abnormal situation, allergic reaction or discomfort occurs during use.
- This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.
- This product is not suitable to use on body parts that are embedded with silicone, metal or other foreign objects. Do not use directly on the eyes.
- It is recommended to wait at least 6 months after any facial treatments such as a face lift or Botox injections before using this product. If you are unsure, consult with your doctor.
- Thoroughly clean the product after use.
- Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
- This product has a special lacquer - when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.
- Repair, removal or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact an authorized dealer to obtain after-sales service if the product has a problem.
- Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
- Do not use this product in locations with high temperatures or high humidity levels (e.g., bathrooms); ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
- Do not insert foreign objects or liquids into the product; this could damage the product.
- Refrain from dropping the product or bumping it against other objects. Avoid damage, collision and excessive pressure by always handling the product with care.
- Do not place this product near strong magnetic or electric fields - this may cause product failure to occur.
- Do not use this product in hazardous locations.
- The product may explode if disposed of in a fire.
- Do not use internally and avoid using near center of the throat to avoid thyroid stimulation.
- Do not use on open sores, wounds or infected areas
- Never use internally (i.e. in mouth, nose, etc.)
- Do not use if pregnant, have a pacemaker, heart disease history or other medical condition.

### Children

Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

### Disposal of Electrical & Electronic Equipment

This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

### Battery

It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the best possible performance. A rechargeable battery can be charged and discharged over one hundred times. However, it will eventually wear out—this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations without contacting an authorized Lifetrans representative.

### Limited Warranty

If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

1. The original invoice or sales receipt with the date of purchase is presented with the defective unit.
2. Lifetrans does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
3. The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrans authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced.

Lifetrans reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.

4. The warranty does not cover product failures which have been caused by use of accessories or other peripheral devices which are not Lifetrans branded original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental breakage.
5. The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
6. Any defective unit or part shall become property of the manufacturer.

### Additional Information

Register your Lifetrans Beauté product online at [www.lifetrans.com](http://www.lifetrans.com).

Please visit our website [www.lifetrans.com](http://www.lifetrans.com) or contact our customer service center at customer [cs@lifetrans.com](mailto:cs@lifetrans.com) to learn more about our new products and special offers.

Product specifications are subject to change at the sole decision and discretion of Lifetrans Switzerland AG, Lifetrans International Ltd. or Li-tek without prior notice. Any devices or accessories not listed under "What's Inside" are in reference to compatibility only and are not included. Product images are for reference purposes only, color may vary from the actual product.

"Lifetrans Beauté", "Lifetrans Beauté Switzerland", "Lifetrans Beauté Professional Skincare", "Lifetrans Switzerland", "Beauté Switzerland", "MyB" are trademarks or registered trademarks owned by Lifetrans Switzerland AG, Lifetrans International Ltd. or Li-tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2016 by Lifetrans Switzerland AG, Lifetrans International Ltd. and Li-tek. All rights reserved.

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

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