
Ultra Photon Light

Red & Blue Light Skin Treatment

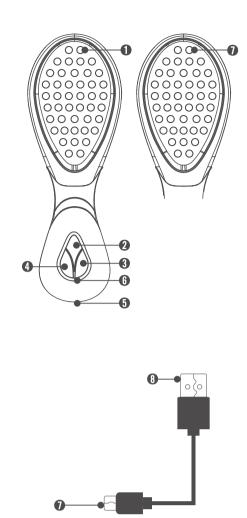EN User Guide

Thank you for purchasing this Lifetrons[®] product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future

How to Use
A. Select one of the two heads
Choose your own color head (1) (red head) or (2) (blue head).

- B. Red light mode

 1. Please use makeup remover and cleanse the skin.

 2. Hold the main unit on the bottom with right hand and pinch the \$\mathbb{0}\$ by the left hand. Load the unit on the top with the head and hear the pop sound when it fits.

 3. Wear the goggles.

 4. Press \$\mathbb{0}\$ blue LED is always on) Di to start, \$\mathbb{0}\$ is on, defaults to the last selected working mode.

 5. Move the \$\mathbb{0}\$ from about 5 cm to the face, and move slowly from bottom to top.

C. Blue light mode

- C. Blue light mode

 1. Please use makeup remover and cleanse the skin.

 2. The right hand holding the host, left hand pinch into the host, pull it out, hear pop sound once been sucked
- place.
 3. Wear the goggles.
 4. Press @ (blue LED is on) with Di sound to start, @ always
- on, the default mode is selected.

 5. Hold away from the facial skin about 5 cm, slowly move it from bottom to top.

- D. Treatment time setting

 1. In the standby mode, press ②, the first level blue LED will be on. After 1 minute, the device will beep three times consecutively.

 2. In the first level position, press ②, the second level blue LED will be on. After 3 minutes, the device will beep for three times consecutively.

 3. In the 2 nd level position, press ②, the third level blue LED will be on. After 5 minutes, the device will beep for three times consecutively.
- three times consecutively.

 4. In the third level position, press ②, the fourth level blue LED will be on. After 12 minutes, the device will beep for
- three times consecutively.

 5. In the fourth position, press ②, the L level blue LED will be on. After 120 minutes, the device will beep for three
- times consecutively.

 6. In the L position, press 2 to enter standby mode. After 1 minute, the device will beep for three times consecutively

E. Pulse light mode

1. In any operating mode, press ③ once (blue LED is on) to enter pulsed mode, ❶ blinks.

2. Press ⑤ (blue LED goes off) again for a "Di" sound, you can return to the working mode, ❶ is on.

F. Shutdown

1. In any operating mode, press , it will beeps three times consecutively and then turn off.

2. In any operating mode, after setting the timing, the device will beep three times consecutively after closing it.

G. Charging 1. Connect () (mini USB plug) and () (mini USB connector). 2. Connect () (USB plug) to a powered USB port (for example, an AC adapter with an USB power supply).

H. Battery Indicator Low Battery: (§ (Green LED) blinks During charging: (§ (red LED) is on Battery full: (§ (red LED) goes out

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EN Safety Precautions

- **Safety Precautions**Please follow the rules outlined in this guide closely to ensure
- Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

 This product is only intended for use on the face.

 Discontinue using the product immediately if any abnormal situation, allergic reaction or discomfort occurs during use.

 This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.

- with your doctor before use.

 Keep the device clean after using.

 Do not expose to eyes directly.

 Before operation, wear eye goggles.

 Do not use any detergent or liquid to clean the surface of this product. Clean the LEDs surface with a dry cloth when needed.

 This product has a special lacquer when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.
- product. Clean the LEDs surface with a ray cloth when needed.

 This product has a special lacquer when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.

 Repair, removal or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact an authorized dealer to obtain after-sales service if the product has a problem.

 Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.

 Do not use this product in locations with high temperatures or high humidify levels (e.g., bathrooms); ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.

 Do not insert foreign objects or liquids into the product; this could damage the product.

 Refrain from dropping the product or bumping it against other objects. Avoid damage, collision and excessive pressure by always handling the product mear strong magnetic or electric fields this may cause the product failure.

 Do not pake this product in hazardous locations.

 The product may explode if disposed in a fire.

 Ould the product with other people.

 Children

ep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

Disposal of Electrical and Electronic Equipment

This symbol indicates that all electrical and electronic
equipment included with this product should not be treated as
household waste. Instead, it should be left at the appropriate
collection point for the recycling of electronic equipment. By
ensuring this product is disposed of correctly, you are helping to
prevent negative consequences for the environment and human
health. To recycle this product, please contact your local city
office, environmental protection agency, household waste
disposal service or where you purchased the product to learn
more about appropriate recycling methods.

Battery
It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the best possible performance.

A rechargeable battery can be charged and discharged over one hundred times. However, it will eventually wear out—this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of it's life span. Please do

not make repairs, modifications, or installations without contacting an authorized Lifetrons representative.

- contacting an authorized Lifetrons representative.

 Limited Warranty
 If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

 1. The original invoice or sales receipt with the date of purchase is presented with the defective unit.

 2. Lifetrons does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.

 3. The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrons authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrons reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.

 4. The warranty does not cover product failures which have been
- caused by use of accessories or other peripheral devices which are not Lifetrons branded original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental

Please visit our website www.lifetrons.com or contact or

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

and discretion of Lifetrons Switzerland As, Interrons internauonal Ltd. or Li-Tek without prior notice. Any devices or accessories no listed under "What's Inside" are in reference to compatibility only and are not included, Product images are for reference purposes only, color may vary from the actual product.

only, color may vary from the actual product.

"Lifetrons Beauté", "Lifetrons Beauté Switzerland", "Lifetrons Beauté Professional Skincare", "Lifetrons Switzerland", "Beauté Switzerland", "MyB" are trademarks or registered trademarks owned by Lifetrons Switzerland AG, Lifetrons International Ltd. or Li-Tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. "Copyright 2017 by Lifetrons Switzerland AG, Lifetrons International Ltd. and Li-Tek. All rights reserved.

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customer service center at customer **cs@lifetrons.com** to learn more about our new products and special offers.

Product specifications are subject to change at the sole decision and discretion of Lifetrons Switzerland AG. Lifetrons International

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- 4. The warranty does not cover product failures which have been
- breakage.

 5. The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.

 6. Any defective unit or part shall become property of the manufacturer.
- Additional Information
 Register your Lifetrons Beauté product online at www.lifetrons.com.

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