





Ultrasonic Facial

With Ions and Thermo Massage



A

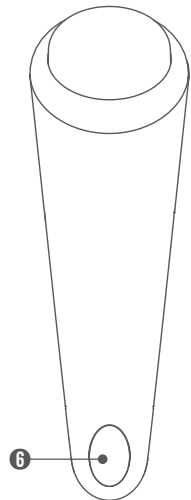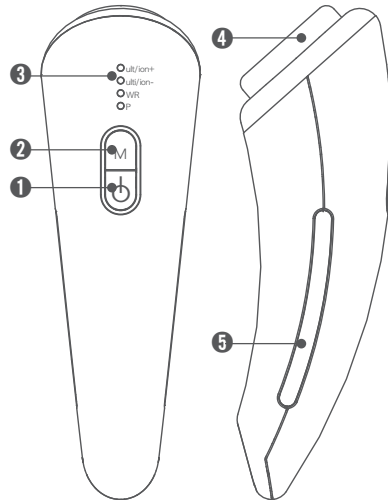

EN User Guide

Thank you for purchasing this Lifetrans® product. This user guide outlines the operating procedures and safety precautions. Please read it carefully to make sure you are using your product properly and retain it for future reference.

How to Use

A. Ultrasonic & Positive Ion Cleansing Mode (ult/ion+)

1. Use on clean skin and be sure to remove makeup before use.
2. Press the power button ❶ to begin. The Ult/ion+ LED ❸ will light up and a beep will sound.
3. Hold the device while touching both conductive side poles ❺ with your hand to activate the ions. The treatment will start and you can place the device head ❹ directly onto your skin to begin. Gently glide the device over wet skin to export impurities from your pores using the power of positive ions which attract dirt. After use, clean your skin again with a toner or similar product to make sure any unseen impurities are cleaned away after treatment.

B. Ultrasonic & Negative Ion Absorption Mode (ult/ion-)

1. Use on clean skin and be sure to remove makeup before use. Serum, moisturizer, cream or facial sheet mask (not included) must be applied in combination with the steps below to ensure proper operation and maximum efficacy.
2. Press the power button ❶ to begin. Press the M button ❷ once to switch to the Ult/ion- treatment mode. The Ult/ion- LED ❸ will light up and once beep will sound.
3. Hold the device while touching both conductive side poles ❺ with your hand to activate the ions. The treatment will start and you can place the device head ❹ directly onto your skin to begin. Gently glide the device over your skincare or mask on your face to help boost the absorption of nutrients and enhance your skincare efficacy. Also may help to improve skin's texture and reduce signs of aging such as fine lines.

C. Warming Micro-Vibration Massage Treatment (WR)

1. Use on clean skin and be sure to remove makeup before use. Serum, moisturizer or cream (not included) must be applied in combination with the steps below to ensure proper operation and maximum efficacy.
2. Press the power button ❶ to begin. Press the M button ❷ twice to switch to the WR treatment mode. The WR LED ❸ will light up and two beeps will sound.
3. The treatment will start and you can place the device head ❹ directly onto your skin to begin. Gently glide the device over your skincare or mask on your face for a soothing, relaxing and warming massage. This may help relieve fatigue, relax facial muscles and increase blood circulation.

D. How to Care for the Device

Thoroughly clean the product with a dry cloth after use. Do not immerse in water or use harsh cleaning chemicals.

E. How to Shut Down the Device

The device will automatically shut down after 10 minutes of continuous use. If you want to shorten the usage time, simply press and hold the power button ❶ for 3 seconds. The device will beep 3 times before shutting down.

F. How to Recharge the Device

1. Open the charging port cap ❹ and connect the micro-USB cable to the charging port.
2. Connect the other side of USB cable with a power supply (5V e.g., Computer USB port or mobile charger).

G. LED Indicator

Low battery: P LED ❸ flashes every 5 seconds (only when the device is in use)

Charging: P LED ❸ flashes every 2 seconds

Fully charged: P LED ❸ light is no longer visible

EN Safety Precautions

Safety Precautions

Please follow the rules outlined in this guide closely to ensure your safety and flawless operation:

- This product is only intended for use on the face.
- Discontinue using the product immediately if any abnormal situation, allergic reaction or discomfort occurs during use.
- This product is a skin care device and is not made for medical purposes. If you are under medical treatment, please consult with your doctor before use.
- This product is not suitable to use on body parts that are embedded with silicone, metal or other foreign objects. Do not use directly on the eyes.
- It is recommended to wait at least 6 months after any facial treatments such as a face lift or botox injections before using this product. If you are unsure, consult with your doctor.
- Thoroughly clean the product after use.
- Do not use any kind of liquid to clean the product. Keep in mind that the surface of the product should only be wiped with a clean dry cloth.
- This product has a special lacquer - when cleaning, do not use any corrosive solution which may cause damage to the lacquer or appearance.
- Repair, removal or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact an authorized dealer to obtain after-sales service if the product has a problem.
- Be cautious of any edges, uneven surfaces, metal parts, accessories or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
- Do not use this product in locations with high temperatures or high humidity levels (e.g., bathrooms); ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
- Do not insert foreign objects or liquids into the product; this could damage the product.
- Refrain from dropping the product or bumping it against other objects. Avoid damage, collision and excessive pressure by always handling the product with care.
- Do not place this product near strong magnetic or electric fields - this may cause the product failure.
- Do not use this product in hazardous locations.
- The product may explode if disposed in a fire.

Children

Keep out of reach of children. Do not allow children to play with your product or its accessories. They could hurt themselves or others and could accidentally damage the product. Your product and its accessories may contain small parts that could become detached and create a choking hazard.

Disposal of Electrical and Electronic Equipment

⚠ This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

Battery

It is recommended that you completely use the initial charge up, then recharge the product at least two more times to ensure the best possible performance.

A rechargeable battery can be charged and discharged over

one hundred times. However, it will eventually wear out—this is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations without contacting an authorized Lifetrans representative.

Limited Warranty

If within 12 months of purchase this product or any part thereof is proven to be defective by reason of faulty workmanship or defect in materials, we will repair or replace the same product free of charge of labor or materials on condition that:

1. The original invoice or sales receipt with the date of purchase is presented with the defective unit.
2. Lifetrans does not take any responsibility for lost information of any kind and will not reimburse you for any such loss.
3. The limited warranty is only valid if the original proof of purchase for this product, issued by a Lifetrans authorized dealer, specifying the date of purchase, is presented with the product to be repaired or replaced. Lifetrans reserves the right to refuse warranty service if this information has been removed or changed after the original purchase of the product from the dealer.
4. The warranty does not cover product failures which have been caused by use of accessories or other peripheral devices which are not Lifetrans branded original accessories intended for use with the product. The warranty does not cover any damage caused by misuse, abuse, neglect or accidental breakage.
5. The warranty does not cover product failures caused by modifications, repairs or disassembly by any person not authorized by the manufacturer.
6. Any defective unit or part shall become property of the manufacturer.

Additional Information

Register your Lifetrans Beauté product online at www.lifetrans.com.

Please visit our website www.lifetrans.com or contact our customer service center at customer cs@lifetrans.com to learn more about our new products and special offers.

This user guide was originally composed in English and was subsequently translated into other languages. The fidelity of subsequent translations cannot be guaranteed. In case of conflict between the English version and another language version, the English version takes precedence.

Product specifications are subject to change at the sole decision and discretion of Lifetrans Switzerland AG, Lifetrans International Ltd. or Li-Tek without prior notice. Any devices or accessories not listed under "What's Inside" are in reference to compatibility only and are not included. Product images are for reference purposes only, color may vary from the actual product.

"Lifetrans Beauté", "Lifetrans Beauté Switzerland", "Lifetrans Beauté Professional Skincare", "Lifetrans Switzerland", "Beauté Switzerland", "MyB" are trademarks or registered trademarks owned by Lifetrans Switzerland AG, Lifetrans International Ltd. or Li-Tek. All other brands and product names listed herein are trademarks or registered trademarks owned by their respective holders and are hereby acknowledged and recognized as such. ©Copyright 2018 by Lifetrans Switzerland AG, Lifetrans International Ltd. and Li-Tek. All rights reserved.